

# sinus- SuperFlex-DS

SOLUTIONS FOR THE HYPOPLASTIC  
LEFT HEART SYNDROME (HLHS)



CONGENITAL HEART DISEASES

# THE HYPOPLASTIC LEFT HEART SYNDROME




sinus-SuperFlex-DS



Congenital Heart Disease (CHD) is defined as a structural abnormality of the heart and/or great vessels that is present at birth. The prevalence of CHD is around 0.8 – 1 % of all live births.<sup>1</sup> CHD accounts for almost 28 % of deaths from birth defects in the first month of life and approximately 50 % of the deaths during the first 12 months.<sup>2</sup>



# HYPOPLASTIC LEFT HEART SYNDROME

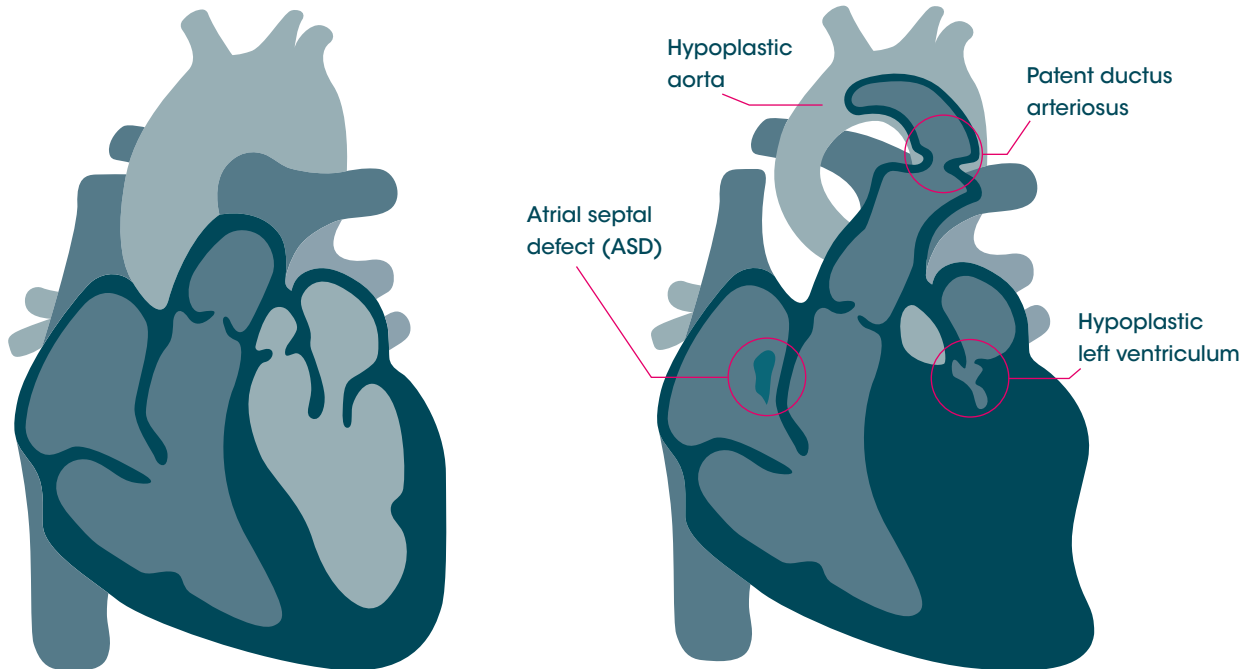

Normal heart

Hypoplastic Left Heart Syndrome

Schematic drawing 1: Comparison of a normal heart anatomy versus Hypoplastic Left Heart Syndrome


Hypoplastic Left Heart Syndrome (HLHS) is the term used for a variety of simultaneously occurring malformations of the heart and aorta in neonates, that are classified within the category LVOTO (Left ventricular out-flow tract obstruction) together with the CHD subtypes coarctation of the aorta and aortic stenosis<sup>3</sup>.

HLHS has a prevalence of 0.178 per 1000 birth and represents 2.56% of all CHD subtypes.<sup>4</sup>

HLHS is characterized by an underdeveloped left ventricle that is completely closed off from the aorta, so that blood flow to the aorta is provided by the ductus arteriosus to ensure that essential organs are supplied with blood (see schematic drawing 1).

HLHS is among the most serious congenital heart defects and is fatal in 99% of cases if left untreated.

# THERE ARE TWO DISTINCT APPROACHES TO TREAT HLHS:

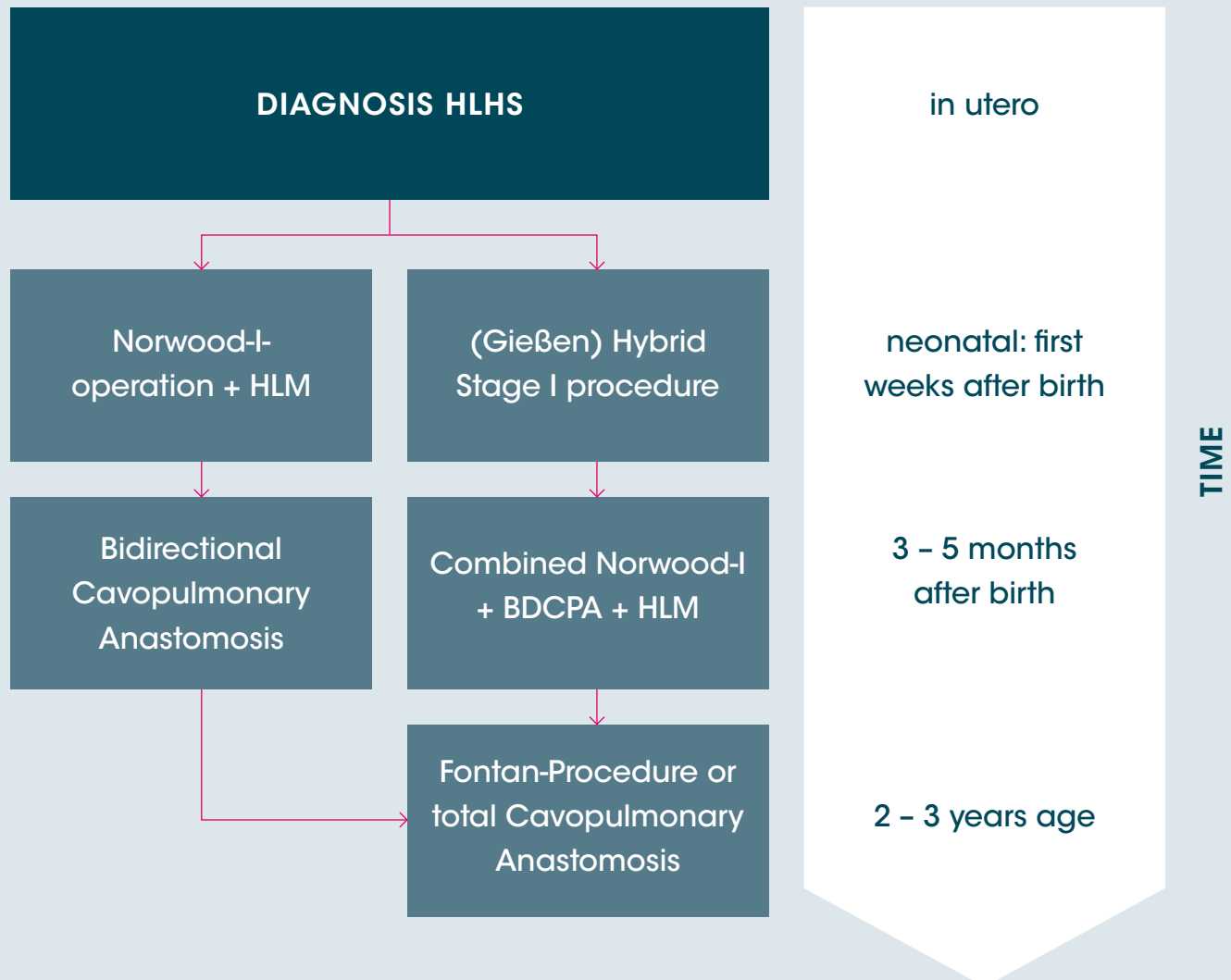

Sources: 1 M, Skovránek J (1989) Prevalence, treatment, and outcome of heart disease in live-born children: a prospective analysis of 91,823 live-born children. *Pediatr Cardiol* 10:205-211)

2 (Benjamin EJ, et al; American Heart Association Statistics Committee and Stroke Statistics Subcommittee. Heart Disease and Stroke Statistics-2017 Update: A Report From the American Heart Association. *Circulation*. 2017 Mar 7;135(10):e146-e603. doi: 10.1161/CIR.0000000000000485. Epub 2017 Jan 25. Erratum in: *Circulation*. 2017 Mar 7;135(10 ):e646. Erratum in: *Circulation*. 2017 Sep 5;136(10 ):e196. PMID: 28122885; PMCID: PMC5408160.)

3 Classification of CHD subtypes based on Bottos's method (2007).

4 (Yingjuan Liu, Sen Chen, Liesl Zühlke, Graeme C Black, Mun-kit Choy, Ningxiu Li, Bernard D Keavney, Global birth prevalence of congenital heart defects 1970-2017: updated systematic review and meta-analysis of 260 studies, *International Journal of Epidemiology*, Volume 48, Issue 2, April 2019, Pages 455-463, <https://doi.org/10.1093/ije/dyz009>)


5



# THE HYBRID PROCEDURE

The hybrid procedure is a cooperative procedure combining techniques of catheter intervention and cardiac surgery.

sinus-SuperFlex-DS

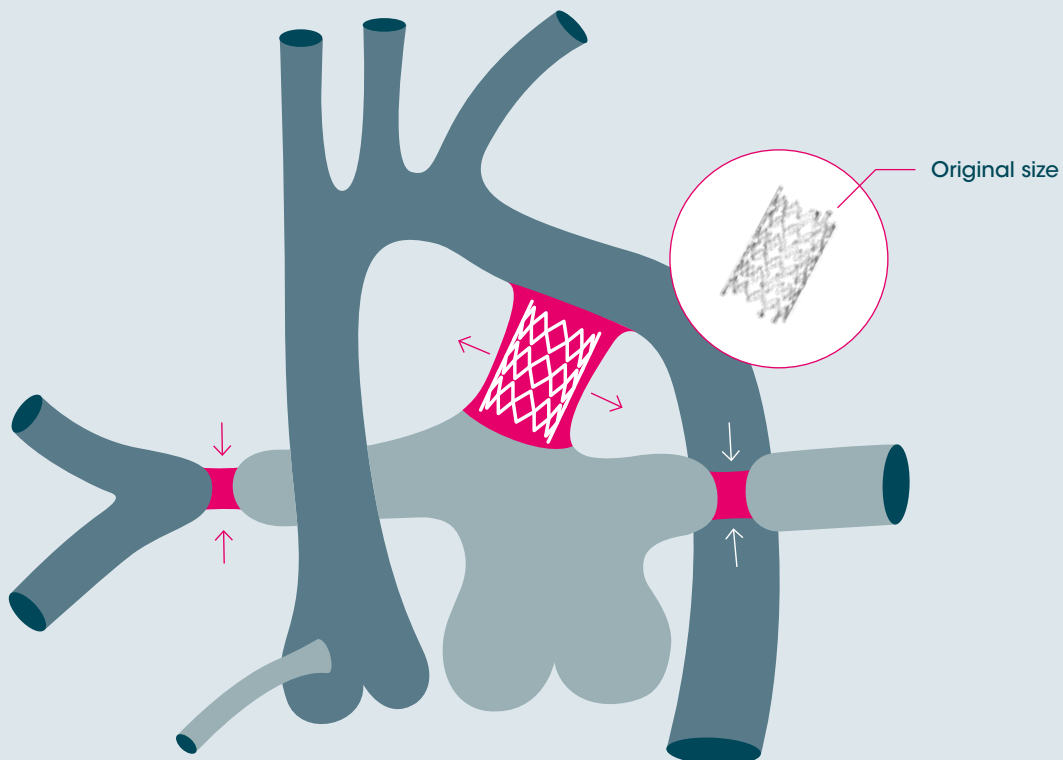

After birth, both branches of the pulmonary artery are surgically narrowed by reducing the high flow of blood to the lungs and directing the blood through the patent ductus arteriosus into the aorta. A stent is implanted to maintain blood flow through the ductus arteriosus.

The open-cell design of the sinus-SuperFlex-DS stent is designed to provide the flexibility and radial force required to maintain patency of the ductus arteriosus until its surgical removal.

## ADVANTAGES OF THE HYBRID INTERVENTION FOR THE NEONATAL / INFANT PATIENT – AT A GLANCE

---

**The Hybrid Intervention avoids a surgical intervention involving a Heart-Lung-Machine (HLM) very early after birth:**

The complex reconstruction of the aorta does not take place immediately after birth but only at an age of 4–6 months – when the child is more stable.

---

**The Hybrid approach is overall regarded as less invasive to the neonatal / infant patient:**

There is NO complete circulatory or cardiac arrest required thus reducing the risk for neurological damage and cardiac muscle damage.



---

A close-up photograph of a hand holding a baby's foot, symbolizing care and support. The hand is positioned in the foreground, with the fingers gently cradling the foot. The background is blurred, showing more of the hand and the foot. The overall tone is warm and intimate.

**EXCEPTIONAL  
SOLUTIONS  
FOR EXCEPTIONAL  
SITUATIONS.**



The self-expanding sinus-SuperFlex-DS vascular stent system is specifically designed to provide temporary patency of the ductus arteriosus. It is indicated to improve or maintain duct dependent systemic circulation in neonates and infants with hypoplastic left heart syndrome/complex (HLHS/HLHC).

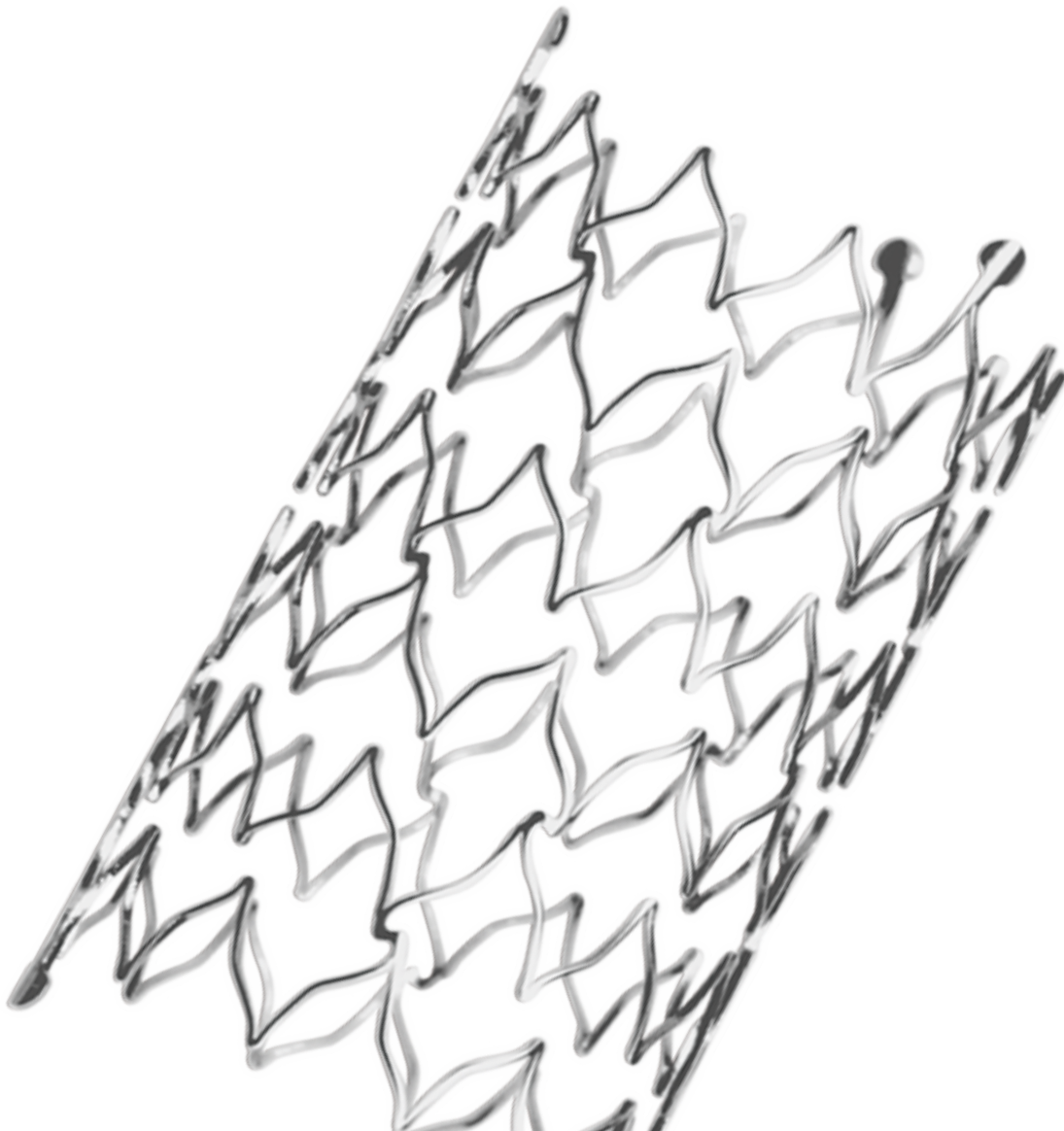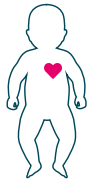

# UNIQUE DESIGN FEATURES

## sinus-SuperFlex-DS

- The **first and only dedicated Ductus Arteriosus stent** on the market
- Available in an **extensive product range** which allows in the vast majority of patient morphologies for the entire duct to be covered with one single stent
- Open-cell system for a **very high level of flexibility** and an optimal vessel adaptation
- Tantalum markers designed for a good radiopacity and **precise stent placement**
- Electro-polishing process for particular smooth surfaces and **reduced thrombogenicity**

The **sinus-SuperFlex-DS** stent has been specifically designed for the treatment of **Hypoplastic Left Heart Syndrome** and has **Class III approval (93/42/EEC)**. The stent design is based on the proven sinus technology, in which the stents are **manufactured seamlessly** from one piece of nitinol, completely without junctions or soldered joints. This approach makes them particularly **resistant to stent breakage** and gives them a **very high radial force even with the shortest stent lengths**. A electropolishing process ensures particularly smooth surfaces and additionally **reduces thrombogenicity**<sup>1</sup>.

## OPTIMAL PATIENT SAFETY

A 4F covered delivery system ensures a simple and gentle access route via the femoral artery. Thus no need for a long vascular sheath or a vascular sheath traumatising of the cardiac structure<sup>2</sup>

Lower risk of stent embolization due to the open-cell stent design<sup>2</sup>

Stent release is possible without balloon dilatation

No risk of complete interruption of the blood flow in the duct which is often fatal<sup>2</sup>

No risk of embolisation caused by retraction of poorly deflated balloons when using balloon expandable stents<sup>2</sup>

Sources: 1 <https://pubmed.ncbi.nlm.nih.gov/16095686/>

2 Mienert, Tine: Das hypoplastische linke Herz: Gießen-Hybrid-Prozedur Stage I mit Einsatz eines speziell entwickelten, selbstexpandierenden Ductus-Stents

Arterial



maximum  
.018 inch  
guide wire



4F Ø  
delivery system



Nitinol  
stent



Braided  
sheath



Radiopaque  
markers



Anti-Jump  
Technique™



Open-cell  
design



# MAXIMUM FLEXIBILITY

The highly flexible 4F application system of the sinus-DS stents provides two access options:



**TRANSVENOUS:**  
over the femoral vein




**TRANSARTERIAL:**  
over the femoral artery

## DETAILS sinus-SuperFlex-DS

### LENGTHS




### DIAMETERS

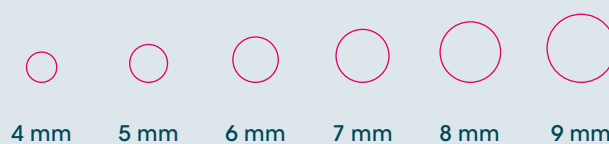


The application system features an exceptionally easy-to-use pull-back-system. Radiopaque markers, braided sheath and atraumatic soft tip facilitate guidance and positioning. To ensure accurate positioning, the stent is fixed in the application device by the anti-jump system – an integrated mechanism which prevents the stent from jumping.



Atraumatic  
soft tip

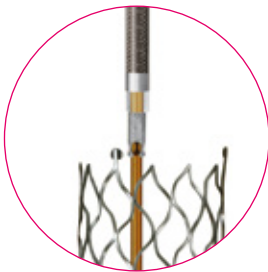

Anti-Jump  
Technique™

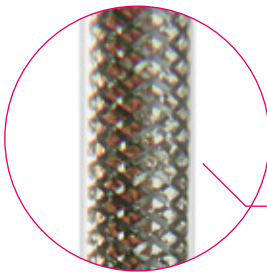

Braided  
outer sheath

4 – 9 mm Ø

## ORDER CODES

	Lengths					
	12	15	18	20	22	24
Ø4	–	8804-2015	8804-2018	8804-2020	8804-2022	8804-2024
Ø5	–	8805-2015	–	8805-2020	–	8805-2024
Ø6	–	8806-2015	–	8806-2020	–	8806-2024
Ø7	8807-2012	8807-2015	8807-2018	8807-2020	8807-2022	8807-2024
Ø8	8808-2012	8808-2015	8808-2018	8808-2020	8808-2022	8808-2024
Ø9	–	–	–	8809-2020	8809-2022	8809-2024



4F / 85 cm  
delivery  
system



Adapted to  
0.018 inch  
guide wire




Tantal marker  
distal 4 /  
proximal 4



Box /1 unit




# DUCTAL STENTING PROCEDURE

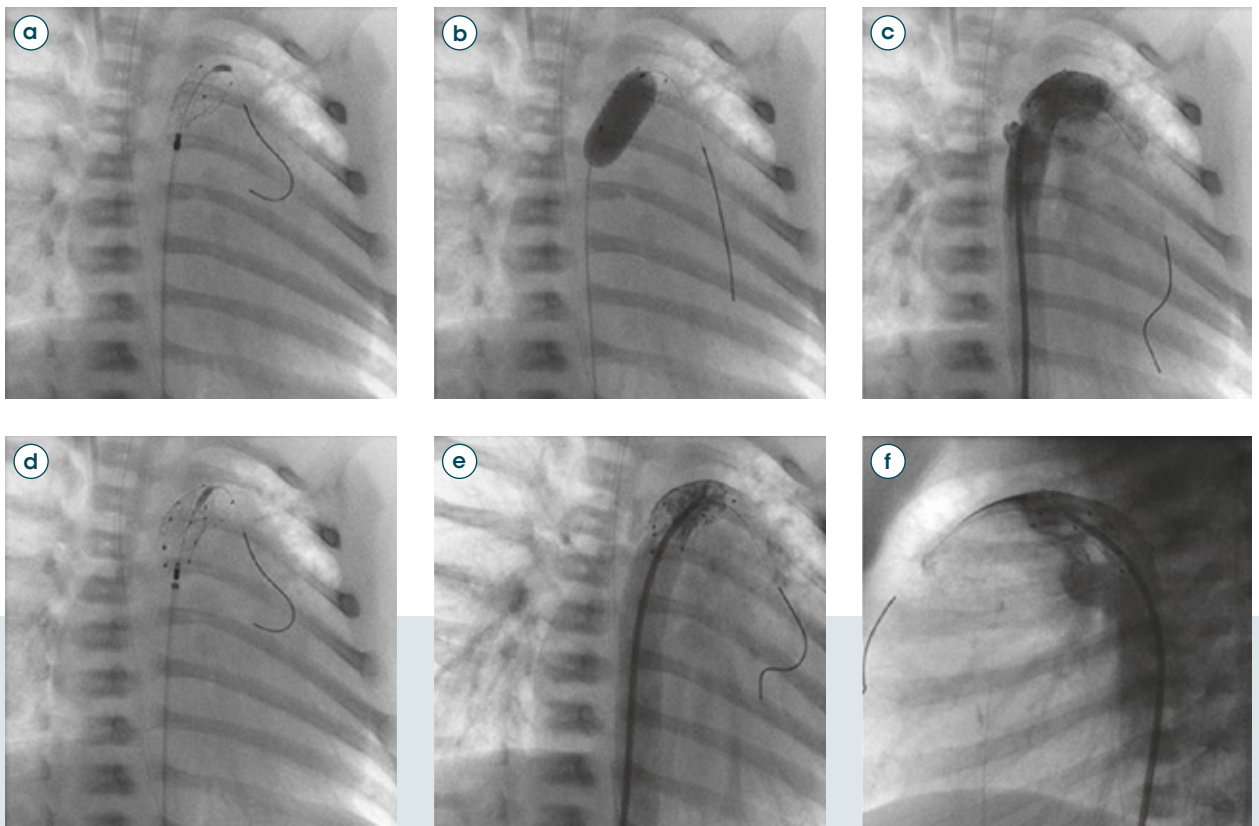

- (a) The sinus-SuperFlex-DS stent (8x18mm) shows incomplete expansion at its pulmonary end after placement into the Ductus Arteriosus.
- (b) Subsequent dilation of the stent using a pediatric valvuloplasty catheter (8x20mm) results in complete expansion of the stent.
- (c) Angiographic assessment reveals inadequate coverage of the Ductus Arteriosus.
- (d) Placement of a second sinus-SuperFlex-DS stent (8x15mm) implanted in an overlapping technique.
- (e) (f) Favorable angiographic result complete stent expansion and total coverage of the Ductus Arteriosus with sufficient safety margins

The sinus-SuperFlex-DS stent before implantation and seven weeks after implant removal – the stent was easily removed from the patient not requiring any dissection and leaving smooth tissue behind.

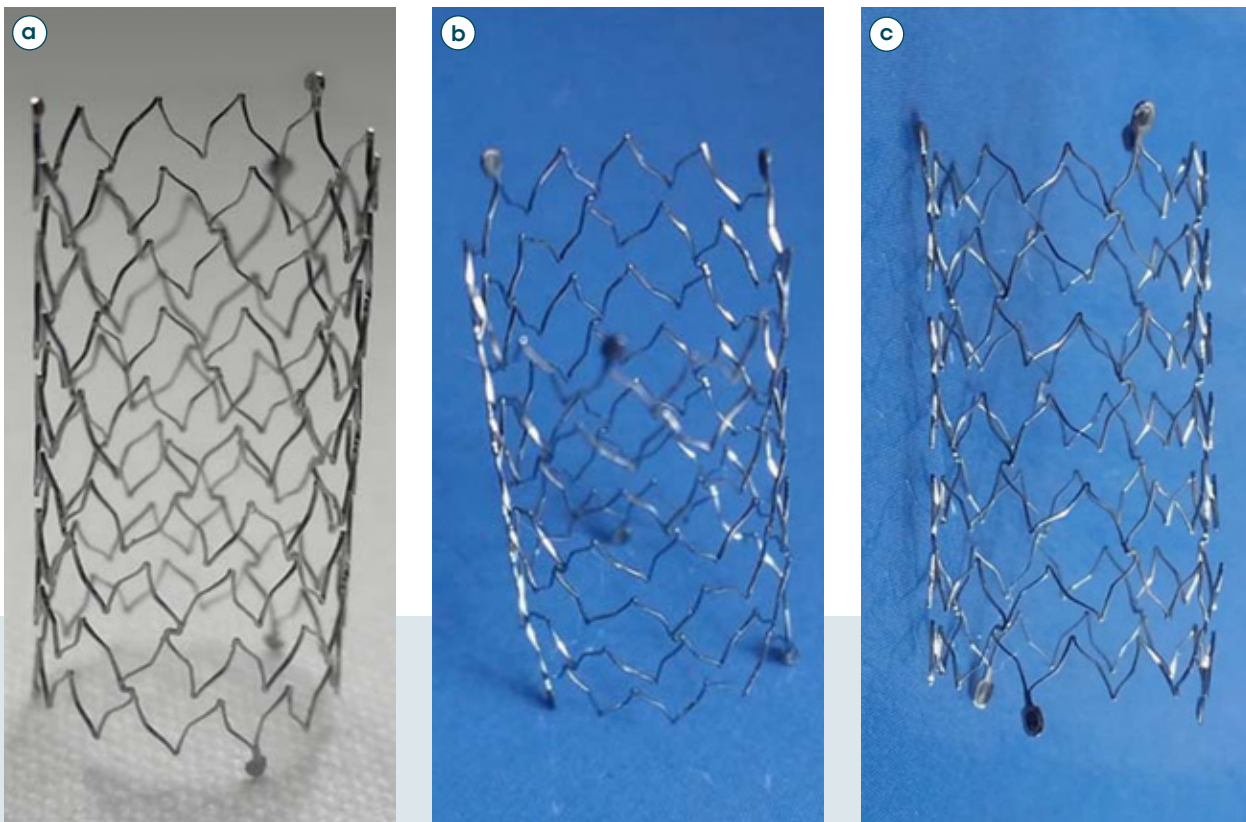

- Ⓐ An example of a self-expanding, open-cell, nitinol sinus-SuperFlex-DS stent.
- Ⓑ and Ⓒ 7 weeks after hybrid ductal stenting, the stent was easily slid out not requiring any dissection and leaving smooth tissue behind.

Sources:

Left page: Butera, Gianfranco & Chessa, Massimo & Eicken, Andreas & Thomson, John. (2019). Atlas of Cardiac Catheterization for Congenital Heart Disease. 10.1007/978-3-319-72443-0.

Chapter 29. Dietmar Schranz and Hakan Akintuerk: Hypoplastic Left Heart Syndrome: The Giessen Hybrid Approach. Page 255 - 261.

This page: Goreczny, Sebastian & Qureshi, Shakeel & Rosenthal, Eric & Krasemann, Thomas & Nassar, Mohamed & Anderson, David & Morgan, Gareth. (2016). Comparison of self-expandable and balloon-expanding stents for hybrid ductal stenting in hypoplastic left heart complex. Cardiology in the Young. 27. 1-9. 10.1017/S1047951116001347.

# YOUR PATIENT IS THE FOCUS OF OUR EVERYDAY ACTIONS

Optimed, based in Ettlingen (Germany), has been engaged in research, development, manufacturing and worldwide distribution of high-end medical products for minimal invasive therapy since 1996. More than 200 employees in our headquarters, production facilities and our global sales teams ensure every day that we provide you and your patients with the highest quality medical devices.



Visit our website at [www.optimed.com](http://www.optimed.com) to discover our comprehensive portfolio and to view our latest product catalogue.

